



# CR0q

Children Bed

SK-A2



**SAIKANG** / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE

Stock Code: 870098

## TECHNICAL PARAMETERS

External L*W	1600*740*1040mm
Bed Height	500mm
Mattress Platform	1580*660mm

## TECHNICAL CONFIGURATION

5" Covered Castors	4pcs
Full Length Side Rail	1set
Intergrated Bed Ends	1set
IV Pole	1pc
IV Pole Prevision	4pcs
Drainage Hook	2pcs

## STANDARD ACCESSORY

One telescopic IV pole

## OPTIONAL ACCESSORY

Basic foam mattress - 12cm thickness  
Anti-decubitus foam mattress - 14cm thickness  
Vertical oxygen holder  
Single arm lifting pole  
Utility Shelf

## QUALITY CONTROL



## RICH PARTS

LEADING BRAND OF  
CHINA MEDICAL FURNITURE



Punching platform, smooth surface no welding scar, bearing 250KG load. Option resin X-ray translucent platform



Powder coated steel frame, MDF material fire-proof plate. size of the head&foot board can be customized. With carton photograph




Foldable guardrail with protection handle, modern and beautiful, easy to install, easy to clean, anti-pinch, protect patient's safety



TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts  
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500time



Full SS IV pole, bearing 15kg

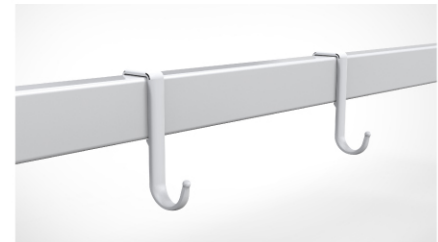

Plastic drainage hook, movable  
Shoes holder: 800\*400mm, adjustable  
Oxygen cylinder rack - adjustable for anysize tank  
Lifting pole - Bearing 765N safety





## INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

### JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901   
[www.chinasaikang.com](http://www.chinasaikang.com)

#### SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,  
Jiangsu Province,China.  
TEL: +86-512-58525850  
FAX: +86-512-58520713  
E-mail: [inquiry@jssaikang.com](mailto:inquiry@jssaikang.com)

#### SAIKANG (USA Branch)

Saikang Medical Technology, Inc  
17800 Castleton St Ste 180,  
City of Industry, Ca 91748  
Tel: +1 949648-2776

#### SAIKANG (UK Branch)

Uk Healthward International Limited  
Owner(s) Address DEPT 114,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +44 2032399738

#### SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.  
Viraj Complex Unit 3 Mombasa Road,  
PO Box 14262-00800, Nairobi.  
TEL: +254 733-704 690

#### SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC  
X22-S03-England Cluster,  
International City, Dubai UAE  
TEL: +20 128 8962497



Contact



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.  
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.