



# CX2x

## TWO FUNCTIONS CHILDREN BED

Registration code: SK-C2




**SAIKANG** Medical

Healthier  
your life

# CX2x

## FEATURES AND HIGHLIGHTS

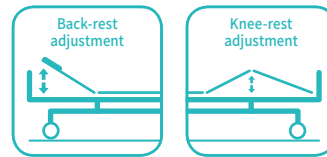

### Technical parameters

▪ External size (LxWxH)	1880x900x1260 mm
▪ Bed height	585 mm
▪ Mattress platform	1770x720 mm
▪ Back-rest adjustment	0-75° (±10°)
▪ Knee-rest adjustment	0-35° (±10°)



### Technical configuration

<input checked="" type="checkbox"/> Manual Cranks	2pcs
<input checked="" type="checkbox"/> 5" Double Side Castors	4pcs
<input checked="" type="checkbox"/> Full Length Side Rail	1set
<input checked="" type="checkbox"/> Plastic Bed Ends with Cartoon Bear	1pc
<input checked="" type="checkbox"/> IV Pole	1pc
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs
<input checked="" type="checkbox"/> Buffer Air Spring	4pcs



**Design**

The bed frame is formed by sheet metal bending and welding, which reduces the dead angle of the plastic spraying and makes the product more beautiful.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding**

Panasonic robotic ensure, 360° full smooth wedding.

**Hand Control**

The Bed With Two Function, Controlled By manual cranks.

**Platform**

The concave vent hole design of the plastic spray bed panel has strong air permeability and can effectively prevent bedsores. The bed panel is connected by sheet metal connectors, which is more convenient for assembly..

**Guardrail**

Adopt aluminum alloy liftable guardrail, the distance between the protective tubes is less than 8 cm, and there are 3 gears to adjust the height, anti-pinch hands to protect children's safety.

**Head & Foot Board**

The bed ends are integrally formed by PP blow molding, which is warm and beautiful and can effectively protect the bed and the patient from colliding with the wall.

**Castors**

5" silent casters, each with independent brakes.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV Pole**

Full SS IV pole, bearing 15kg.

**Shoes Holder**

Utility holder & Shoes holder.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue