



SKB039(D)

AMBULANCE STRETCHER TROLLEY

Registration code: SKC-3B5



SAIKANG / Medical

Healthier
your life

SKB039(D)

FEATURES AND HIGHLIGHTS



Technical parameters

▪ Full size (LxWxH)	1900x550x920mm
▪ Lowest size (LxWxH)	1900x550x250mm
▪ Back-rest lifting	0-75°
▪ Net weight	35kg
▪ Load bearing	≤159kg

Safety Lock

The safety lock is locked when transporting a patient, and unlocked when getting on an ambulance.



Caster

5" silicone caster. Advanced protection universal brake wheel, easy to move, simple to operate, reliable performance.



IV Pole

Telescopic stainless steel infusion rod, strong bearing capacity, more convenient to use.



Function

The mechanical folding legs are controlled by the handles at both ends of the stretcher, which can be easily lowered when entering the ambulance.



Back Rest Lifting

Back is adjustable which makes the wounded feel comfortable. Back rest lifting 0-75°.



Fixing Devices

Can be locked by the fixing devices when on the ambulance car.




Structure

Adopt high-quality aluminum alloy material, fitting aluminum alloy die-casting, light weight, anti-corrosion and anti-rust.



Cushion

Adopt sponge cushion to make the injured and sick lying comfortably.

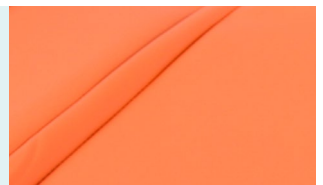

Leather / Safety Belts

Artificial leather, bright appearance, soft texture, abrasion resistance, folding resistance, acid and alkali resistance. With 2 safety belts, make patient safer.



Guardrail

The aluminum alloy flip guardrail is convenient for patients to get on and off the stretcher, and is equipped with two safety belts.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue