




SKB041-2

Patient
Transportation
Trolley



SAIKANG / Medical

LEADING BRAND OF CHINA MEDICAL FURNITURE
Stock Code: 870098

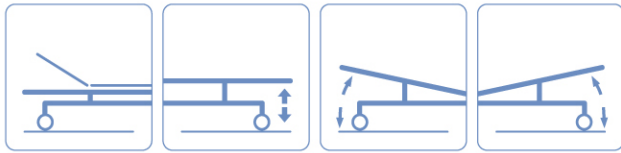

Back-rest Height Adjustment Reverse Trendelenburg Trendelenburg

MAIN TECHNICAL PARAMETERS

| | |
|-----------------------|----------------------|
| Size | 2120*800*(630-930)mm |
| Back-rest | 0-75° |
| Reverse Trendelenburg | 0-15° |
| Trendelenburg | 0-15° |

TECHNICAL CONFIGURATION

| | |
|---------------------------------|-------|
| Hydraulic Pump | 2pcs |
| Back Rest Air Spring | 1pc |
| ABS Side Rail | 1set |
| 8" Castors with Central Locking | 4pcs |
| Directional Wheel | 1pc |
| IV Pole Previsions | 2pcs |
| IV Pole | 1set |
| Central Locking Break Pedal | 2sets |
| Transfer Mattress | 1pc |
| Directional Pedal | 2pcs |

QUALITY CONTROL

METAL MATERIAL

Baosteel from Fortune Global 500
Panasonic Robotic ensure 360° full
smooth wedding

500 | 360°

PLASTIC

Environmental plastic can
be sterilized in 100°C,
tensile upto 30MPa

100°C | 30MPa

STEELS | PAINTING

11 Process epoxy painting,
ASTM testing anti-bacterial, paint
thickness 0.12mm, brightness
60°, paint can resist 50kg impact

0.12mm | 60° | 50kg

BRAKE


Break spring keep qualified
after 10,000times
uninterupted break

HYDRAULIC PUMPS

Height adjustment,
trendelenburg and
revised trendelenburg
through two USA
hydraulic pumps

CENTER WHEEL

Can be engaged to facilitate
direction control, enabling
users to navigate the
stretcher easily in hectic and
dynamic environments.

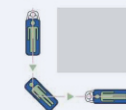

RICH PARTS

LEADING BRAND OF
CHINA MEDICAL FURNITURE



Hydraulic treadle for rise



- ① Leg section lows down to make head higher (Than Leg Section)
- ② Back section lows down to make leg section higher (Than head Section)
- ③ Step on the middle of the pedal, make the bed base fall down steadily



Four double side casters of 6 inch diameter; simultaneously lockable by the step pedal at both sides.

TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts



Full SS IV pole, bearing 15kg



Water proof mattrass with 2 belts





INNOVATION OF SAIKANG

LEADING BRAND OF CHINA MEDICAL FURNITURE



Design Innovation



Research Innovation



Manufacturing Innovation



Maintenance Innovation



Service Innovation

JIANGSU SAIKANG MEDICAL EQUIPMENT CO.,LTD.

24hrs: 0086-18021231901 
www.chinasaikang.com

SAIKANG (China Headquarter)

No.35 Lehong Road,Zhangjiagang City,
Jiangsu Province,China.
TEL: +86-512-58525850
FAX: +86-512-58520713
E-mail: inquiry@jssaikang.com

SAIKANG (USA Branch)

Saikang Medical Technology, Inc
17800 Castleton St Ste 180,
City of Industry, Ca 91748
Tel: +1 949648-2776

SAIKANG (UK Branch)

Uk Healthward International Limited
Owner(s) Address DEPT 114,
196 High Road Wood Green,
London, N22 8HH
United Kingdom
TEL: +44 2032399738

SAIKANG (East Africa Branch)

Saikang Medical Solutions (East Africa) Ltd.
Viraj Complex Unit 3 Mombasa Road,
PO Box 14262-00800, Nairobi.
TEL: +254 733-704 690

SAIKANG (Mid East Branch)

Saikang Medical Supplies (Mid East) LLC
X22-S03-England Cluster,
International City, Dubai UAE
TEL: +20 128 8962497



Catalogue



Website

As the circumstance of printing, the photos might have minor difference on color, tactile from real products.
As the update of design and improvement, parameters and outlook of products might change, thanks for understanding.