



SKE-180

ELECTRIC
DIALYSIS CHAIR



SAIKANG / Medical


Healthier
your life

SKE-180

ELECTRIC DIALYSIS CHAIR

Technical parameters

▪ Size	1800×580mm
▪ Back	770×580mm
▪ Seat	470×580mm
▪ Leg	380×540mm
▪ Foot board	410×200mm
▪ Seat height	550-750mm
▪ Armrest	500×160mm
▪ Mattress thickness	60mm
▪ Back rest lifting	-8°- 70°
▪ Leg rest lifting	20° -75°
▪ Armrest swing out	0°-270°

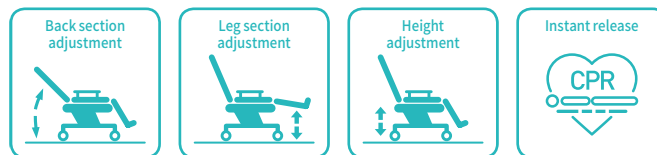

- Voltage and frequency: 100-240V · 50/60 Hz
- Plug as required (hospital grade)

Technical configuration

<input checked="" type="checkbox"/> With CPR	
<input checked="" type="checkbox"/> Hand controller	1pc
<input checked="" type="checkbox"/> Power cord	1pc
<input checked="" type="checkbox"/> 3" silent caster	1set
<input checked="" type="checkbox"/> Small dining table	1pc



Design

European design, 4 sections ergonomics design.

**Steels**

Baosteel from Fortune Global 500.

**Plastic**

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

**Welding**


Panasonic Robotic ensure 360° full smooth wedding.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Armrest**

Armrest updated, can be universal adjust. User can adjust any position he/she wants. Armrest is R120° concave curve, avoid donor's arm sliding, ensure the user be safety and comfortable.

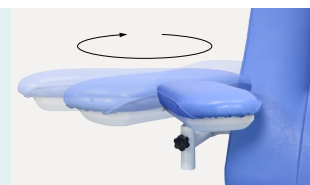**Leather**

Optional genuine leather or supebric cover,with artisanal making,comfortable hand feeling.

**Foot Board**

Foot pedal connects with leg rest, adjustable, can automatic extension the chair lengthh 250mm, realize both foot step and leg support functions.

**Table Board**

Adjustable ABS tray table, with antiskid convex edge, tank for cup and cellphone. USB/5V optional. Very convenient for donor.

**IV Pole**

Full stainless steel IV pole, bearing 15 kg.



ULTRAFIBER LEATHER COLOR TABLE



Off white 5509-1



Off white 5509-2



Beige 5509-3



Pink 5510-45




Light rose red 2892-15



Rose red 2201-30



Wine red 2209-22

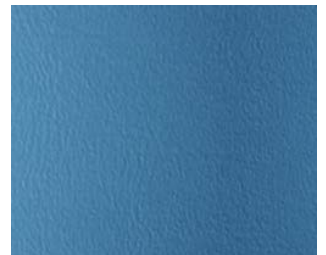

Light blue 5510-43



Sky blue 5510-7



Dark lake blue 2201-40



Light grey 5509-14



Gray 5510-12



Lake blue 5509-27



Navy blue 2201-42



Light green 2210-88



Green 5510-20

Note: The above is the physical scan , there may be a certain color difference in kind , are to prevail in kind.



CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue