



SKP006

AIR MATTRESS



SKP006

FEATURES AND HIGHLIGHTS

Mattress

▪ Arrangement	2.5" Bubble pad
▪ Material	Medical Grade PVC
▪ Dimension	2000×900×70mm
▪ Weight	3kg

Pump

<input checked="" type="checkbox"/> Electrical	AC220-240V/100-120V,50HZ/60HZ
<input checked="" type="checkbox"/> Air Output	4L/min
<input checked="" type="checkbox"/> Pressure Range	12kPa(90mmHg)
<input checked="" type="checkbox"/> Cycle Time:	5-6min



Design

The bubble alternately fluctuate, alternately inflate and deflate, cycling every 6 minutes to prevent bedsores without turning over.




Bubble Mattress

The elastic airbags are closely arranged, and the stress points on the body are more uniform and dispersed. It is ventilated and breathable, dry and comfortable, and the air circulates naturally, not stuffy, to avoid the growth of bacteria.



Material

Medical material, soft and comfortable, breathable and skin-friendly, waterproof and dirt-resistant, environmentally friendly and tasteless.

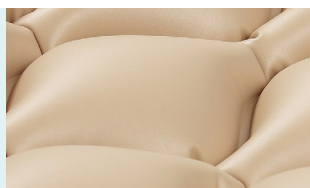

Smart Air Pump

Silent design is less than or equal to 45 decibels, intelligent pressure, you can freely adjust the softness and hardness of the air cushion according to the needs.



Angle

It can be bent freely from multiple angles to meet the needs of turning over. It can also be folded away for storage.



Airway

Double airway circulation, wave massage is more comfortable, circulation ventilation, natural air circulation.



Tightness

Good sealing, integrated heat-sealing, no cracking or air leakage.



Energy saving

Stable and reliable, energy-saving and power-saving, stable operation for 24 hours, only one kilowatt hour for 5 days.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue