



Grace Care series

# D8d

## MULTIFUNCTIONAL ELECTRIC BED

Registration code: SKD-C




**SAIKANG** / Medical

Healthier  
your life

# D8d

## FEATURES AND HIGHLIGHTS

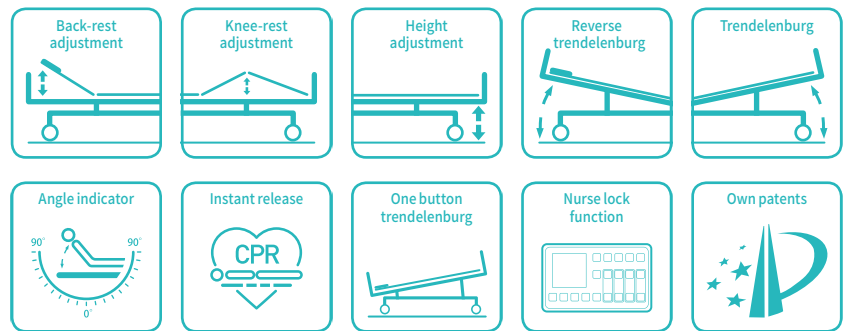

### Technical parameters

▪ External size (LxW):	2280x1090 mm
▪ Mattress platform:	1925x900 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	435-725 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Reverse Trendelenburg:	0-12°
▪ Trendelenburg:	0-12°
▪ Angle indicator:	0-90°
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	



### Technical configuration

<input checked="" type="checkbox"/> Manual CPR	
<input checked="" type="checkbox"/> One button stop	
<input checked="" type="checkbox"/> One button trendelenburg	
<input checked="" type="checkbox"/> Nurse lock system	
<input checked="" type="checkbox"/> Angle indicator for back-rest and Trendelenburg position	4 pcs.
<input checked="" type="checkbox"/> Linear motor	4 pcs.
<input checked="" type="checkbox"/> 5" double side castors with central locking system	4 pcs.
<input checked="" type="checkbox"/> ABS side rail with integrated control panel	1 set
<input checked="" type="checkbox"/> ABS bed ends with integrated LCD control panel	1 set
<input checked="" type="checkbox"/> φ85 mm wall bumper wheel	4 pcs.
<input checked="" type="checkbox"/> Trapeze socket	2 pcs.
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs.
<input checked="" type="checkbox"/> Drainage hook	2 pcs.
<input checked="" type="checkbox"/> Central locking pedal	2 pcs.

### Optional

<input type="checkbox"/> Electrical CPR	
<input type="checkbox"/> Weighing system with LCD Screen	
<input type="checkbox"/> Bed extension	
<input type="checkbox"/> Back-up battery	
<input type="checkbox"/> Cardiac position	
<input type="checkbox"/> ABS bed surface	
<input type="checkbox"/> X-ray bed surface	
<input type="checkbox"/> Monkey pole	
<input type="checkbox"/> Balkan traction frame	
<input type="checkbox"/> Directional wheel	
<input type="checkbox"/> Basic foam mattress (12 cm thickness)	
<input type="checkbox"/> Anti-decubitus foam mattress (14 cm thickness)	
<input type="checkbox"/> Vertical oxygen holder	
<input type="checkbox"/> Single arm lifting pole	
<input type="checkbox"/> Utility shelf	

**Head and Foot Board**

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Bed Surface**

Stamping forming platform, smooth surface no welding scar, bearing 250 kg load. Option ABS platform and resin X-ray translucent platform.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding. Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Steels | Painting**

Baosteel from Fortune Global 500. 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Virgin ABS**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

**LCD**

Unique design head foot board with built-in controller. 4.3 inch LCD display, shows moving angle. Button can endure 100,000 times press.

**Side Rail**

Tuck away side rail with buffer, built-in controller and angle indicator.

**Wheels**


TPR tire no worn out after running 30 km, save anti-winding hard shell, united forming without bolts. Pass dynamic test: Bearing 120 kg, run 30 km and pass obstacles 500 time.

**Durable Brakes**

Break spring keep qualified after 10,000 times uninteruped break.

**Angle Indicator**

Patented angle indicator.

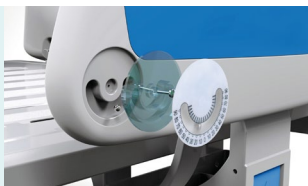**Bumper Wheels**

Patented bumper wheel detachable, max. protect the bed and patient during transportation.

**Safety**

IEC60601 Safty Standard.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic hook.

**Optional**

Battery: 18 times movement with 100 kg bearing, 130 times movement without bearing.  
Nurse control board.  
Storage rack under foot board drawer.  
Store controll board or linens.  
Utility basket: big size 600x540x80 mm.







## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### USA

**Saikang Medical Technology, Inc.**  
17800 Castleton St Ste 180,  
City of Industry, 91748, California.  
Tel.: +1 949648 2776

### United Kingdom

**UK Healthward International Ltd.**  
Dept. 111, 196 High Road Wood  
Green, N22 8HH, London.  
Tel.: +44 2032399738

### Russia

**Saikang Medical Russia LLC**  
Truda 174,  
454080 Chelyabinsk.  
Tel.: +7 (351) 214 5566 / 8811 / 1516

### East Africa

**Saikang Medical Solutions Ltd.**  
Viraj Complex Unit 3, Mombasa  
Road, P.O. Box 14262-00800, Nairobi.  
Kenya.  
Tel.: +254 733 704 690

### Middle East

**Saikang Medical Supplies LLC**  
X22-S03 England Cluster,  
International City, Dubai.  
United Arab Emirates.  
Tel.: +20 128 8962497

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue