



Oasis Series

Y6w

THREE FUNCTIONS ELECTRIC BED

Registration code: SKD-C

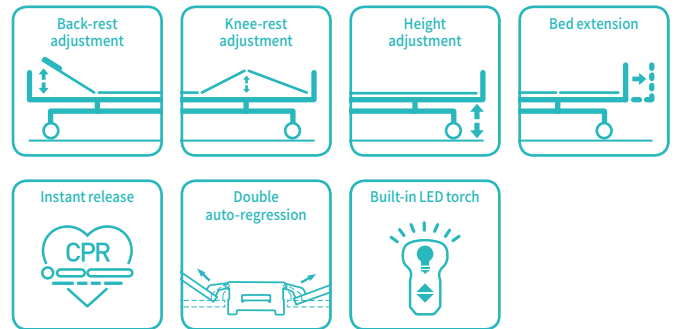


SAIKANG / Medical

Healthier
your life

Y6w

FEATURES AND HIGHLIGHTS



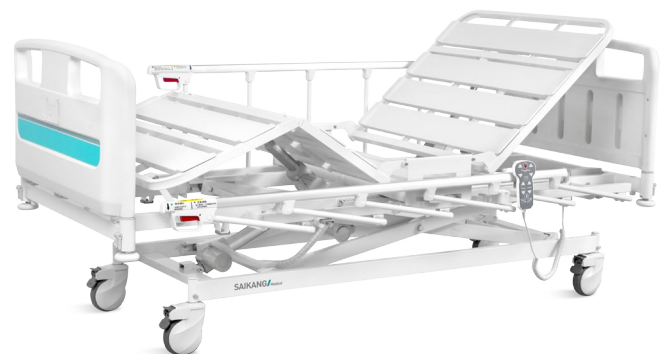
Technical parameters

▪ External size (LxW):	2130x1015 mm
▪ Mattress platform:	1925x880 mm
▪ Safe working load:	250 kg
▪ Height adjustment:	365-805 mm
▪ Back-rest adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-35° (±10°)
▪ Bed Extension:	200mm
▪ Auto-regression:	100+60mm
▪ Voltage and frequency:	100-240 V · 50/60 Hz
▪ Plug as required (hospital grade)	



Technical configuration

<input checked="" type="checkbox"/> Linear Motor	3pcs
<input checked="" type="checkbox"/> 5" Hardshell Castors	4pcs
<input checked="" type="checkbox"/> Foldable Side Rail with ProErgo-Handle™	1set
<input checked="" type="checkbox"/> ABS Oasis Bed Ends with Safe Lock	1set
<input checked="" type="checkbox"/> PP Bed Platform, Removable & Washable	9pcs
<input checked="" type="checkbox"/> Ø100mm Wall Bumper Wheel	4pcs
<input checked="" type="checkbox"/> IV Pole	1set
<input checked="" type="checkbox"/> IV Pole Prevision	4pcs
<input checked="" type="checkbox"/> Drainage Hook	2pcs
<input checked="" type="checkbox"/> ABS Mattress Stopper	2pcs
<input checked="" type="checkbox"/> Hand Controller	1pc



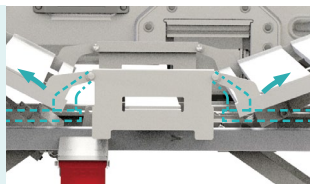
Optional

<input type="checkbox"/> Manual CPR	
<input type="checkbox"/> Electric CPR	
<input type="checkbox"/> Bed Extension 20cm	
<input type="checkbox"/> Back-up Battery	
<input type="checkbox"/> Nurse Control Panel	
<input type="checkbox"/> Emergency Stop Button	
<input type="checkbox"/> Directional Wheels Function	
<input type="checkbox"/> Monkey Pole	
<input type="checkbox"/> Balkan frame	



Double Regression Function

Auto-Regression 10+6 cm
When repositioning the backrest and thigh rest, the Y8y makes the back plate retreat 10cm and the thigh plate retreats 6cm so to minimise the pressure and shear forces in the patient's pelvic region. This results in a significant reduction in the development of pressure ulcers.

**Electric Motor**

Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500
11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.

**Welding | Chassis**

Panasonic Robotic ensure, 360° full smooth wedding.
Solid bed base made of 3.0mm stamped rotating plate.

**Head & Foot Board**

Patented design Detachable head & foot board, using the eccentric locking structure, one step quick removable.

**Safety Guardrail**

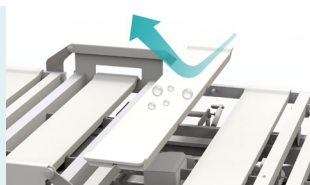
L-shaped foldable guardrail(Patent), bended tube, solid installing, 0 shaking, anti-pinch, protect patient's safety.

**Castors**

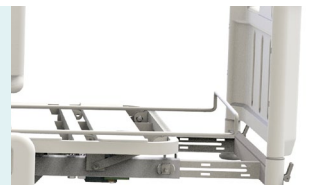
5" hardshell casters
TPR tire no worn out after running 30KM, save anti-winding hard shell, united forming without bolts.
Pass Dynamic test: Bearing 120kg run 30KM and pass obstacles 500 time.

**Joint-Free Technology**

Removable plastic covers formed from 9 removable blow moulded lamella to minimise the slots where unwanted infections may hide.

**Bed Extension**

Extended bed platform up to 20CM, comfortable for the tallest of patients (option).

**Mattress**

Mattress Holder.

**Wall Bumper Wheel**

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.

**IV pole**

Stainless Steel Column with plastic hook.

**Drainage Hooks**

Mobile plastic drainage hooks.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 📞 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue