

Grace Care series

D2w


TWO FUNCTIONS MANUAL BED

Registration code: SK-CD1



D2w

FEATURES AND HIGHLIGHTS

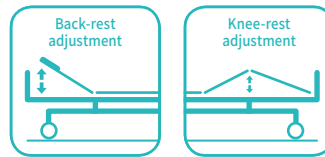

Technical parameters

| | |
|-------------------------|------------------|
| ▪ External size (LxWxH) | 2245x1000x510 mm |
| ▪ Mattress platform | 1925x830 mm |
| ▪ Safe working load | 250 kg |
| ▪ Back-rest adjustment | 0-75° (±10°) |
| ▪ Knee-rest adjustment | 0-35° (±10°) |



Technical configuration

| | |
|---|------|
| <input checked="" type="checkbox"/> Manual Cranks | 2pcs |
| <input checked="" type="checkbox"/> 5" Hardshell Castors | 4pcs |
| <input checked="" type="checkbox"/> Foldable Side Rail with ProErgo-Handle™ | 1set |
| <input checked="" type="checkbox"/> ABS Grace Bed Ends with Safe Lock | 1set |
| <input checked="" type="checkbox"/> IV Pole | 1set |
| <input checked="" type="checkbox"/> IV Pole Prevision | 4pcs |
| <input checked="" type="checkbox"/> Drainage Hook | 2pcs |
| <input checked="" type="checkbox"/> Shoes Holder | 1set |



Optional

| | |
|--|--|
| <input type="checkbox"/> Monkey Pole | |
| <input type="checkbox"/> Balkan Traction Frame | |
| <input type="checkbox"/> Directional Wheel | |



Design

Extra long and height head foot board, with curve surrounded design, provide comforting feeling for patients.

**Plastic**

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa
Virgin ABS material.

**Steels | Painting**

Baosteel from Fortune Global 500 11
Process epoxy painting, ASTM testing
anti-bacterial, paint thickness 0.12
mm, brightness 60°, paint can resist
50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full
smooth welding Solid bed base
made of 5.0 mm stamped rotating
plate.

**Hand Control**

Retractable hand crank, 100,000
times remaining intact
Extreme protection device, sounds
when tipped to the max position,
Strong pallets to increase the bed
height adjustment.

**Platform**

Stamping forming platform, smooth
surface no welding scar, bearing
250KG load. Option resin X-ray
translucent platform.

**Guardrail**

L-shaped foldable guardrail (Patent)
with ProErgo-Handle™, patented parts.
with red and green mark to indicated
if the siderail is locked or not. Bended
tube, solid installing, 0 shaking, anti-
pinch, protect patient's safety.

**Head & Foot Board**

Head and foot board can be reversable,
extra height design, size: 1000*563mm;
ABS Grace Bed Ends with Safe Lock,
automatic lock without any manual
operation on lock, on/off removable
mark on hook lock.

**Castors**

5" hardshell casters
TPR tire no worn out after running
30KM, save anti-winding hard shell,
united forming without bolts.
Pass Dynamic test: Bearing 120kg run
30KM and pass obstacles 500 time.

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic
hook.

**Shoes Holder**

Shoes holder: 800*400mm,
adjustable





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎️ +86 18021231901 · 🕒 24/7

✉️ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue