



CT8k

ELECTRIC CHILDREN BED

Registration code: SKD-C

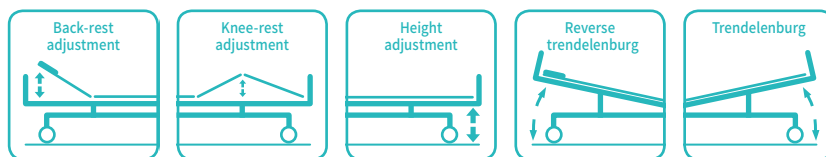


SAIKANG / Medical

Healthier
your life

CT8k

FEATURES AND HIGHLIGHTS



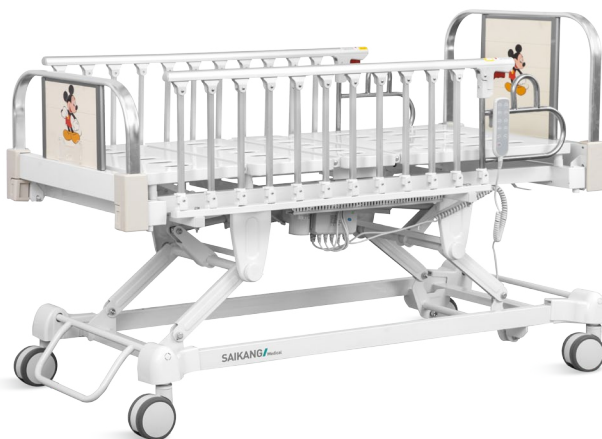
Technical parameters

External size (LxWxH)	1800x780x480-770 mm
Back-rest adjustment	0-75° (±10°)
Knee-rest adjustment	0-35° (±10°)
Reverse Trendelenburg:	0-12°
Trendelenburg:	0-12°
Voltage and frequency:	100-240 V · 50/60 Hz
Plug as required (hospital grade)	



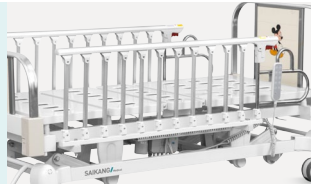
Technical configuration

<input checked="" type="checkbox"/> Linear Motor	4pcs
<input checked="" type="checkbox"/> 5" Double side castors, central locking	4pcs
<input checked="" type="checkbox"/> Aluminum alloy side rail with extra protection	1set
<input checked="" type="checkbox"/> Aluminum alloy bed ends	1set
<input checked="" type="checkbox"/> IV pole	1pc
<input checked="" type="checkbox"/> IV pole prevision	4pcs
<input checked="" type="checkbox"/> Drainage hook	2pcs



Design

Aluminum alloy side rail, 12-staff, providing full protection, meets safety standards and provide comforting feeling for patients.

**Steels | Painting**

Baosteel from Fortune Global 500 11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12 mm, brightness 60°, paint can resist 50 kg impact.

**Welding | Chassis**

Panasonic robotic ensure, 360° full smooth wedding Solid bed base made of 5.0 mm stamped rotating plate.

**Electric Motor**

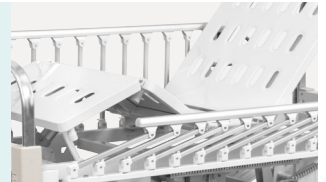
Linear motor, with UL/ROHS/EN standard. Noise under 50 dB, with lifetime 20,000 times movement.

**Head & Foot Board**

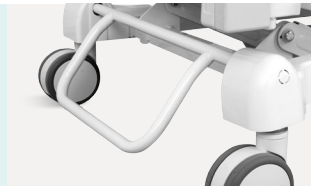
The bed ends are integrally formed by PP blow molding, which is warm and beautiful and can effectively protect the bed and the patient from colliding with the wall.

**Bed Platform**

Stamping forming platform, smooth surface no welding scar, bearing 250KG load.

**Braking System**

Break spring keep qualified after 10,000times uninteruped break.

**Wheels**

5-inch double-sided casters, smooth and smooth running, durable with central locking. Castors with TPR tire, no worn-out after running 30KM, united forming without bolts.

**Bumper Cover**

Reduce the impact of the collision, safe and slinky.

**Nurse Control**

Hook type nurse control.

**IV pole**

Stainless Steel Column with plastic hook.

**Drainage Hooks**

Mobile plastic drainage hooks.





CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



WORLDWIDE BRANCHES

USA

Saikang Medical Technology, Inc.
17800 Castleton St Ste 180,
City of Industry, 91748, California.
Tel.: +1 949648 2776

United Kingdom

UK Healthward International Ltd.
Dept. 111, 196 High Road Wood
Green, N22 8HH, London.
Tel.: +44 2032399738

Russia

Saikang Medical Russia LLC
Truda 174,
454080 Chelyabinsk.
Tel.: +7 (351) 214 5566 / 8811 / 1516

East Africa

Saikang Medical Solutions Ltd.
Viraj Complex Unit 3, Mombasa
Road, P.O. Box 14262-00800, Nairobi.
Kenya.
Tel.: +254 733 704 690

Middle East

Saikang Medical Supplies LLC
X22-S03 England Cluster,
International City, Dubai.
United Arab Emirates.
Tel.: +20 128 8962497

SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ sales@chinasaikang.com

🌐 www.saikangmedical.com



Company video



General catalogue