



# V6k

## ELECTRIC BED

Registration code: SKD-C



**SAIKANG** / Medical

Healthier  
your life

# V6k

## FEATURES AND HIGHLIGHTS



### Technical parameters

External size (LxW):	2085x995 mm
Mattress platform:	1925x900 mm
Safe working load:	200 kg
Height adjustment:	420-730mm
Back-rest adjustment:	0-75° (±5°)
Knee-rest adjustment:	0-46° (±5°)


### Optional

- Battery Back-up
- Bed Extension
- Nurse Controller
- Emergency Stop Button
- Auto-regression Function (Back Section)
- Night Light
- Manual CPR
- X-ray Translucent Mattress Platform (Back Section)
- Utility Shelf
- Trapeze Bracket/Pole

### Technical configuration

<input checked="" type="checkbox"/> Linear motor	3 pcs
<input checked="" type="checkbox"/> Hand controller	1 pc
<input checked="" type="checkbox"/> 5" covered castors	4 pcs
<input checked="" type="checkbox"/> Foldable side rail	1 set
<input checked="" type="checkbox"/> PP bed ends with safe lock with $\Phi 65$ bumper wheels	1 set
<input checked="" type="checkbox"/> PP bed platform with integrated mattress retainer	1 set
<input checked="" type="checkbox"/> IV pole 4 hooks	1 set
<input checked="" type="checkbox"/> IV pole prevision	4 pcs
<input checked="" type="checkbox"/> Drainage hook	2 pcs

### Optional color (MOQ100pcs- color changing)

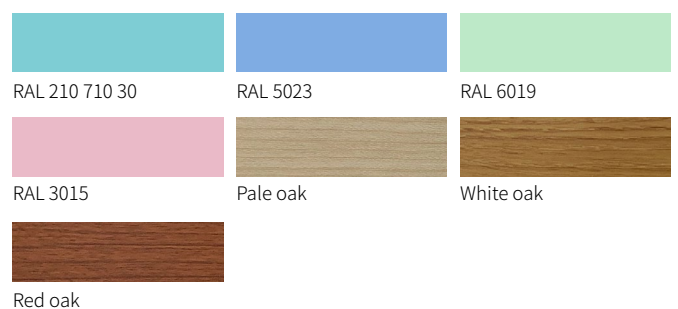

# QUALITY CONTROL

SAIKANG medical, as a pioneer in the field of medical furniture, always choose branded suppliers strictly to ensure the safety and reliability of raw material. SAIKANG adopt robot welding, automatic epoxy coating, efficient assembly line, dedicated operator, rigorous quality control, only to guarantee SAIKANG's commitment on quality



## Design

The unique germ-free design meets hospitals' strict infection control requirement. Seamless PP head and foot board and platform ensure the instant and sufficient sterilization.



## Steels | welding

Baosteel from Fortune Global 500. Panasonic Robotic ensure 360° full smooth welding.



## Plastic | Painting

Environmental plastic can be sterilized in 100°C, tensile upto 30MPa.

11 Process epoxy painting, ASTM testing anti-bacterial, paint thickness 0.12mm, brightness 60°, paint can resist 50kg impact.



## Electric Motor

American motor, with UL/ROHS/EN standard. Noise under 50dB, with lifetime 20,000 times movement.

## Accessories



### Platform

Integrated PP bed platform, adopting advanced manufacturing technology. 4-5 section platform, enables independent installation and easy maintenance.

Enlarge the seat area when back part is raised, providing more comforts for the patients.

Reduce the depression of the sponge at the seat panel.



### Castors

5" covered casters. TPR tire no worn out after running 10 km.

Anti-winding hard shell cover, united forming without bolts.

Pass dynamic test: pass obstacles 40000 times, under the standard 200kg load of the whole bed.

**Head & Foot Board**

LG Chem fresh PP material, Plug in type. Equipped with the bumpers wheel, diameter 65mm, preventing collision.



**Mattress Retainer**

Platform integrated mattress retainer at the both sides of the back and leg section, with a boarder of 33mm height.




**Guardrail**

Aluminum alloy side rail, 6-staff, collapsed type with safe lock.




**Wall bumper wheel**

Horizontal corner bumpers at the outermost periphery of the bed to protect the bed and wall.

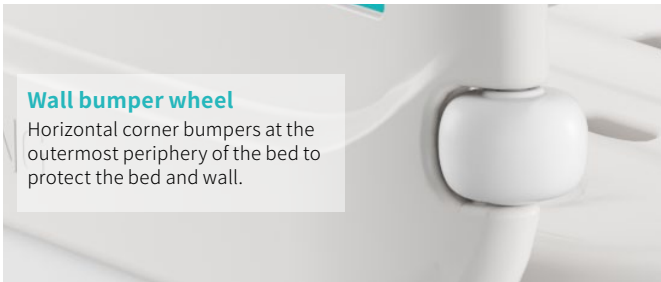

**Drainage Hooks**

Mobile plastic drainage hooks.

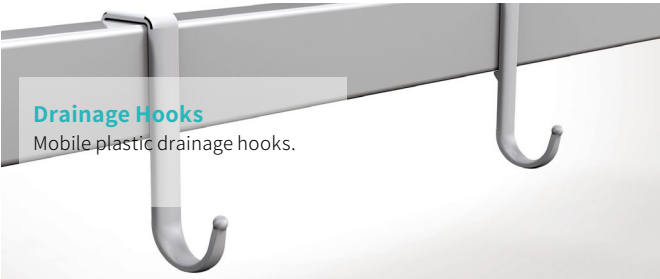


**IV Pole**

SS IV pole, bearing 15kg.



Back-rest adjustment

75°

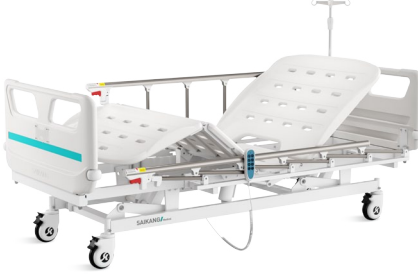

Knee-rest adjustment

46°



Height adjustment

420-730<sub>mm</sub>



## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

**Saikang Medical Technology, Inc**  
W3400 Corrientes, Argentina  
TEL: +5493794312231  
Av. Senador Vergueiro, 4204 - sala 31  
Rudge Ramos - São Bernardo do Campo - SP, Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

**Saikang Medical Technology Ltd**  
Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

**Saikang Medical Technology Ltd**  
Sala Thammassop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

**Uk Healthward International Limited**  
DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

**LLC 'SAIKANG'**  
454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

**Saikang Medical Supplies (Mid East) LLC**  
Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue